



CLEAR BLUE SMILES™

Suresmile® Aligner

#3 Ordering Complete Aligners

Appointment #2

OVERVIEW

Review and confirm for digital setup.

Review ordering protocol for complete clear aligners.

Review clinical protocols for clear aligner treatment.

PREPARATION

Review handout:
#3 Ordering Complete Aligners

Review handout:
#4 Aligner Clinical Protocols

APPOINTMENT

Appointment time commitment: 30 minutes

Contact Ortho Mentor during case triage as needed for additional support.

Appointment #2

STEP 1

Review and confirm modifications in digital setup.

STEP 2

Review Occlusion.

Digital Setup Evaluation:

1. SHO
2. THO

STEP 3

Review ordering complete aligner protocol.

STEP 4

Review Clinical Protocols:

Clear Aligner:

1. Delivery appointment
2. Dispensing protocol
3. Progress visits



CLEAR BLUE SMILES™

Step 1: Review Digital Setup

Step 1: Review Digital Setup

Confirm lab has applied all modifications. Utilize the setup and staging checklist to aid in final treatment plan review. Make any further modifications if needed.

1. Login into SureSmile

Dentsply Sirona

SureSmile®

Log in with your account details

Email address

Password

Login

[Forgot your password?](#)

[Looking to add a new user to your account?](#)

New to SureSmile®?
Click here to create an account

© 2022 OraMetrix. All rights reserved.

[Privacy Policy](#)

2. Select Review setup and staged models

AuraSmile Dental

Patient ID  [New Patient](#)

 Dr. Brian Gaudreault
AuraSmile Dental

 SureSmile Aligner

Sort by: Task Show Tasks of: All Types Tasks for: All View:  

Tasks (4) Patients Shipments Jobs DS CEREC Orders Reports  Active Tasks

Flag	Card	Last name	First name	Patient ID	Owner	Item	Task	Due Date	Notes	Setup Approval Date	Status	Deferred To
<input type="checkbox"/>		Walker	Shekeda	O3P3	Douglas W. Voiers, DDS		Order Retainers  	2022-07-04		2022-06-16	Active	
<input checked="" type="checkbox"/>		Dillon	Sanders	H5Y4	Douglas W. Voiers, DDS		Review setup & staged models 	2022-06-30		2022-07-27	Active	
<input type="checkbox"/>		Wayne	Renee	V9T8	Douglas W. Voiers, DDS		Review setup & staged models 	2022-06-30		2022-07-20	Active	
<input type="checkbox"/>		Miller	Jeff	G6M8	Douglas W. Voiers, DDS		Review setup & staged models 	2022-07-01		2022-07-21	Active	

Setup and Staging Review Opens

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is "Jean Frost (R2N6)" and the current step is "Setup & Staging 1 (Ready for Review)". The main window shows a 3D model of a patient's teeth in a patient view. On the left, a "Checklist and Overview" panel is open, listing several review items, all of which are marked with green checkmarks:

- Setup and Staging Review
- Labial View
- Right Buccal View
- Left Buccal View
- Upper Occlusal
- Lower Occlusal
- IPR / Attachment / Trimline

Below the checklist, there is a "Request Modifications" toggle switch which is currently turned off. At the bottom of the checklist panel, there is a button labeled "Approve & Order Aligners". The bottom of the interface features a navigation bar with playback controls (back, play, forward), a timeline with markers at 1, 5, 10, 15, 20, and 21, and icons for "Legend" and "Treatment Details".

3. Select Overview

The screenshot displays the SureSmile Aligner software interface. At the top, the patient's name "Jean Frost (R2N6)" and the current step "Setup & Staging 1 (Ready for Review)" are visible. A "Checklist and Overview" sidebar on the left contains the following items:

- Overview** (highlighted with a red arrow): Click the Notes Tab and review the notes from the Digital Lab (DL). Go to next step to begin reviewing the setup in different views, and review IPR.
- Setup and Staging Review (checked)
- Labial View (checked)
- Right Buccal View (checked)
- Left Buccal View (checked)
- Upper Occlusal (checked)
- Lower Occlusal (checked)
- IPR / Attachment / Trimline (checked)

Below the checklist is a "Request Modifications" toggle switch. At the bottom of the sidebar is a button labeled "Approve & Order Aligners". The main central area shows a 3D model of a patient's teeth in a "Camera (Patient View) Navigation" perspective. The bottom of the interface features a navigation bar with playback controls (back, play, forward), a progress bar with markers at 1, 5, 10, 15, 20, and 21, and icons for "Legend" and "Treatment Details".

4. Select Notes Tab



Product Notes

All Notes

2022-08-31

05:54:00 **Finished by** Lisette Hitt
OraMetrix #1


2022-08-30

18:29:17 **In Production** Elluany Po...

17:27:18 **Submitted** Dr. Brian G...

17:27:18 **Instructions** Dr. Brian G...

Keep already placed attachments: No
No more than 1mm of cuspid expansion

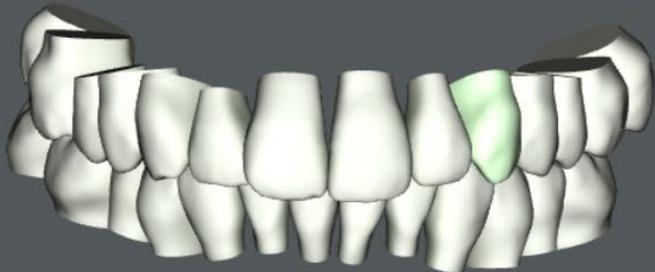

Camera (Patient View) Navigation

Image Assets



Review Digital Lab Notes

The screenshot displays a digital lab interface with three main sections:

- Product Notes:** A sidebar on the left containing a list of notes. A red arrow points to the entry: "05:54:00 Finished by OraMetrix #1 Lisette Hitt". Below this, notes for "2022-08-30" include "18:29:17 In Production", "17:27:18 Submitted", and "17:27:18 Instructions". The instructions note: "Keep already placed attachments: No No more than 1mm of cuspid expansion".
- 3D Model:** A central view showing a 3D digital model of a dental arch. The right side of the model is highlighted in green, corresponding to the "Finished by" note.
- Image Assets:** A panel on the right showing two images: a patient's smiling face and a clinical photograph of the patient's teeth.

At the bottom of the interface, the text "Camera (Patient View) Navigation" is visible.

5. Select Labial View

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is "Jean Frost (R2N6)" and the current stage is "Setup & Staging 1 (Ready for Review)". The main workspace shows a 3D model of the upper dental arch in a labial view. On the left, a "Checklist and Overview" panel lists several views: Overview, Setup and Staging Review, Labial View (highlighted with a red arrow), Right Buccal View, Left Buccal View, Upper Occlusal, Lower Occlusal, and IPR / Attachment / Trimline. Below the checklist is a "Request Modifications" toggle and an "Approve & Order Aligners" button. On the right, an "Image Assets" panel shows a patient's photo and a clinical photo of the teeth. At the bottom, a navigation bar includes playback controls, a timeline from 1 to 21, and buttons for "Legend" and "Treatment Details".

6. Select Right Buccal View

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is Jean Frost (R2N6) and the current stage is Setup & Staging 1 (Ready for Review). The main window shows a 3D model of the teeth with a pink aligner overlay, viewed from the right buccal perspective. A red arrow points to 'Right Buccal View' in the 'Checklist and Overview' panel on the left. The checklist includes: Overview, Setup and Staging Review, Labial View, Right Buccal View, Left Buccal View, Upper Occlusal, Lower Occlusal, and IPR / Attachment / Trimline. To the right, the 'Image Assets' panel shows two clinical photographs of the patient's teeth. The bottom of the interface features a navigation bar with playback controls, a timeline from 1 to 21, and buttons for Legend and Treatment Details.

7. Select Left Buccal View

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is "Jean Frost (R2N6)" and the current stage is "Setup & Staging 1 (Ready for Review)". The main window shows a 3D dental model of the upper and lower teeth with a pink aligner overlay. A red arrow points to the "Left Buccal View" option in the "Checklist and Overview" panel on the left. The checklist includes: Overview, Setup and Staging Review, Labial View, Right Buccal View, Left Buccal View, Upper Occlusal, Lower Occlusal, and IPR / Attachment / Trimline. Below the checklist is a "Request Modifications" toggle and an "Approve & Order Aligners" button. On the right, the "Image Assets" panel shows two clinical photographs of the patient's teeth. At the bottom, there is a navigation bar with a timeline from 1 to 21, a "Legend" button, and a "Treatment Details" button.

8. Select Upper Occlusal View

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is 'Jean Frost (R2N6)' and the current stage is 'Setup & Staging 1 (Ready for Review)'. On the left, a 'Checklist and Overview' panel lists several items, with 'Upper Occlusal' highlighted by a red arrow. The central 3D model shows the upper dental arch with green and yellow highlights on the occlusal surfaces. On the right, an 'Image Assets' panel shows two clinical photographs of the patient's teeth. The bottom of the interface features a navigation bar with a timeline and control icons.

Second Site | Jean Frost (R2N6) | Setup & Staging 1 (Ready for Review) | SureSmile Aligner

Checklist and Overview

- ✓ Overview
- ✓ Setup and Staging Review
 - ✓ Labial View
 - ✓ Right Buccal View
 - ✓ Left Buccal View
 - ▶ Upper Occlusal
 - ✓ Lower Occlusal
 - ✓ IPR / Attachment / Trimline

Request Modifications

Approve & Order Aligners

Image Assets

Camera (Patient View) Navigation

Legend | Treatment Details

9. Select Lower Occlusal View

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is "Jean Frost (R2N6)" and the current stage is "Setup & Staging 1 (Ready for Review)". The main central area shows a 3D model of the lower dental arch with green and yellow highlights on the teeth. To the left, a "Checklist and Overview" panel lists several items, with "Lower Occlusal" highlighted by a red arrow. To the right, an "Image Assets" panel shows two clinical photographs of the lower teeth. At the bottom, there is a navigation bar with a timeline from 1 to 21 and various control icons.

Second Site | Jean Frost (R2N6) | Setup & Staging 1 (Ready for Review) | SureSmile Aligner

Checklist and Overview

- Overview
- Setup and Staging Review
 - Labial View
 - Right Buccal View
 - Left Buccal View
 - Upper Occlusal
 - Lower Occlusal**
 - IPR / Attachment / Trimline

Request Modifications

Approve & Order Aligners

Image Assets

Camera (Patient View) Navigation

Legend | Treatment Details

10. Review IPR/Attachments/Trimline

The screenshot displays the SureSmile Aligner software interface. At the top, the patient's name is "Jean Frost (R2N6)" and the current stage is "Setup & Staging 1 (Ready for Review)". The left sidebar contains a "Checklist and Overview" panel with the following items:

- Overview (checked)
- Setup and Staging Review (checked)
- Labial View (checked)
- Right Buccal View (checked)
- Left Buccal View (checked)
- Upper Occlusal (checked)
- Lower Occlusal (checked)
- IPR / Attachment / Trimline (highlighted with a red arrow)

The main 3D view shows a dental model with a yellow trimline and attachments. A red arrow points to the "IPR / Attachment / Trimline" item in the checklist. The bottom of the interface features a timeline with markers at 1, 5, 10, 15, 20, and 21, and a "Camera (Patient View) Navigation" bar. The right side of the interface shows "Image Assets" with a patient photo and a clinical photo of the teeth.

11. Review Treatment Details

The screenshot displays the SureSmile Aligner software interface for a patient named Jean Frost (R2N6). The main window shows a 3D model of the upper dental arch with a yellow wire and brackets. A checklist on the left indicates that the 'Setup and Staging Review' is complete, including Labial, Right Buccal, Left Buccal, Upper Occlusal, and Lower Occlusal views. The 'IPR / Attachment / Trimeline' view is currently selected. A 'Request Modifications' toggle is turned off. The bottom of the interface features a navigation bar with a timeline from 1 to 21, and buttons for 'Legend' and 'Treatment Details'. A red arrow points to the 'Treatment Details' button.

Second Site | Jean Frost (R2N6) | Setup & Staging 1 (Ready for Review) | SureSmile Aligner

Checklist and Overview

- Overview
- Setup and Staging Review
 - Labial View
 - Right Buccal View
 - Left Buccal View
 - Upper Occlusal
 - Lower Occlusal
 - IPR / Attachment / Trimeline

Request Modifications

Approve & Order Aligners

Image Assets

Camera (Patient View) Navigation

Legend | Treatment Details

Tooth Movements

Stage 29	UR8	UR7	UR6	UR5	UR4	UR3	UR2	UR1	UL1	UL2	UL3	UL4	UL5	UL6	UL7	UL8
mesial (+) / distal (-)		-0.3	-0.1	0.1			-0.1	0.3	-0.3	-0.1	-0.3	-0.2	-0.1	-0.3	-0.3	
buccal (+) / lingual (-)		1.1	-0.3	0.7	0.3	0.5	0.8	-0.3	0.2	0.6	0.2	1.1	0.5			
occlusal (+) / gingival (-)		1.1		-0.1	0.3	0.3	0.4	0.2		0.5	0.2	0.2				
Torque facial (+) / lingual (-)		-16	-4				2		4		-2					
Ang. mesial (+) / distal (-)		6						-4	3	4						
Rot. mesial (+) / distal (-)				-5	5		-7	-18	3	-9	7		-7			
Spacing (+) / Intersection (-)									-0.2							
Do not move	<input type="checkbox"/>															

Cumulative: Cumulative Incremental



Tooth Movements

Stage 29	UR8	UR7	UR6	UR5	UR4	UR3	UR2	UR1	UL1	UL2	UL3	UL4	UL5	UL6	UL7	UL8
mesial (+) / distal (-)																
buccal (+) / lingual (-)																
occlusal (+) / gingival (-)																
Torque facial (+) / lingual (-)																
Ang. mesial (+) / distal (-)																
Rot. mesial (+) / distal (-)																
Spacing (+) / Intersection (-)									-0.2							
Do not move	<input type="checkbox"/>															

Cumulative: Cumulative Incremental




Elastics

Tooth Movements	Elastics	Measurements	IPR Tracking
Class Correction [mm]		Midline Correction [mm]	
0.0 mm		0.0 mm	



Measurements

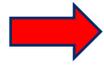

Tooth Movements	Elastics	Measurements	IPR Tracking	
<input checked="" type="radio"/> malocclusion <input type="radio"/> setup				
	Bolton Ratio	Maxilla Sum	Mandible Sum	Lower
All (6..6)	98.1	91.3 mm	89.6 mm	Surplus (6..6): 6.2 mm
Anteriors (3..3)	79.2	44.1 mm	34.9 mm	Surplus (3..3): 0.9 mm

Arch Length Discrepancy		
	Left	Right
Upper	1.0 mm	0.5 mm
Lower	-1.3 mm	-0.3 mm

Measured Overjet
0.0 mm

Measurements

Tooth Movements Elastics **Measurements** IPR Tracking

malocclusion setup 

	Bolton Ratio	Maxilla Sum	Mandible Sum		Lower
All (6..6)	98.1	91.3 mm	89.6 mm	Surplus (6..6):	6.2 mm
Anteriors (3..3)	79.2	44.1 mm	34.9 mm	Surplus (3..3):	0.9 mm

Arch Length Discrepancy		
	Left	Right
Upper	1.4 mm	1.4 mm
Lower	0.4 mm	0.4 mm

Measured Overjet
0.3 mm

IPR Tracking

Tooth Movements		Elastics		Measurements		IPR Tracking										✕
	UR8	UR7	UR6	UR5	UR4	UR3	UR2	UR1	UL1	UL2	UL3	UL4	UL5	UL6	UL7	UL8
Planned								-0.2								
Applied																
Completed	<input type="checkbox"/>															
Completed	<input type="checkbox"/>															
Applied																
Planned				-0.3	-0.3	-0.3	-0.3	-0.3	-0.3							
	LR8	LR7	LR6	LR5	LR4	LR3	LR2	LR1	LL1	LL2	LL3	LL4	LL5	LL6	LL7	LL8

Edit Selection: ▲ ▼ =

Clear Treatment Details Tab

Tooth Movements	Elastics	Measurements	IPR Tracking														
Stage 29		UR8	UR7	UR6	UR5	UR4	UR3	UR2	UR1	UL1	UL2	UL3	UL4	UL5	UL6	UL7	UL8
mesial (+) / distal (-)			-0.3	-0.1	0.1			-0.1	0.3	-0.3	-0.1	-0.3	-0.2	-0.1	-0.3	-0.3	
buccal (+) / lingual (-)			1.1	-0.3	0.7	0.3	0.5	0.8	-0.3	0.2	0.6	0.2	1.1	0.5			
occlusal (+) / gingival (-)			1.1		-0.1	0.3	0.3	0.4	0.2		0.5	0.2	0.2				
Torque facial (+) / lingual (-)			-16	-4				2		4		-2					
Ang. mesial (+) / distal (-)			6						-4	3	4						
Rot. mesial (+) / distal (-)					-5	5		-7	-18	3	-9	7		-7			
Spacing (+) / Intersection (-)										-0.2							
Do not move		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Cumulative:		<input checked="" type="radio"/> Cumulative	<input type="radio"/> Incremental														

Setup and Staging Review Opens

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is "Jean Frost (R2N6)" and the current step is "Setup & Staging 1 (Ready for Review)". The main window shows a 3D model of a patient's teeth in a patient view. On the left, a "Checklist and Overview" panel is open, listing several review items, all of which are marked with green checkmarks:

- Setup and Staging Review
- Labial View
- Right Buccal View
- Left Buccal View
- Upper Occlusal
- Lower Occlusal
- IPR / Attachment / Trimline

Below the checklist, there is a "Request Modifications" toggle switch which is currently turned off. At the bottom of the checklist panel, there is a button labeled "Approve & Order Aligners". The bottom of the interface features a navigation bar with playback controls (back, play, forward), a timeline with markers at 1, 5, 10, 15, 20, and 21, and icons for "Legend" and "Treatment Details".



CLEAR BLUE SMILES™

Step 2: Review SHO and THO

Review SHO

Static hyper-occlusion will appear as red occlusal markings on the final staged models or last set of aligners. This is a planned, intentional collision to obtain a better fit of the occlusion during the treatment planning stage. Inform the patient that it will be adjusted at the aligner delivery appointment.

Select Initial Staged Model


The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is Jean Frost (R2N6) and the current stage is Setup & Staging 1 (Ready for Review). The main window shows a 3D model of a patient's teeth in a patient view. On the left, a 'Checklist and Overview' panel is open, listing various views that have been reviewed, all marked with green checkmarks:

- Overview
- Setup and Staging Review
 - Labial View
 - Right Buccal View
 - Left Buccal View
 - Upper Occlusal
 - Lower Occlusal
 - IPR / Attachment / Trimline

Below the checklist is a 'Request Modifications' toggle switch, which is currently turned off. At the bottom of the interface is a 'Camera (Patient View) Navigation' bar. It includes navigation icons (back, play, forward), a view selector with 'U' and 'L' options, a timeline with a red box around the number '1', and a legend icon. A red arrow points to the '1' on the timeline.

Review SHO


Initial Stage

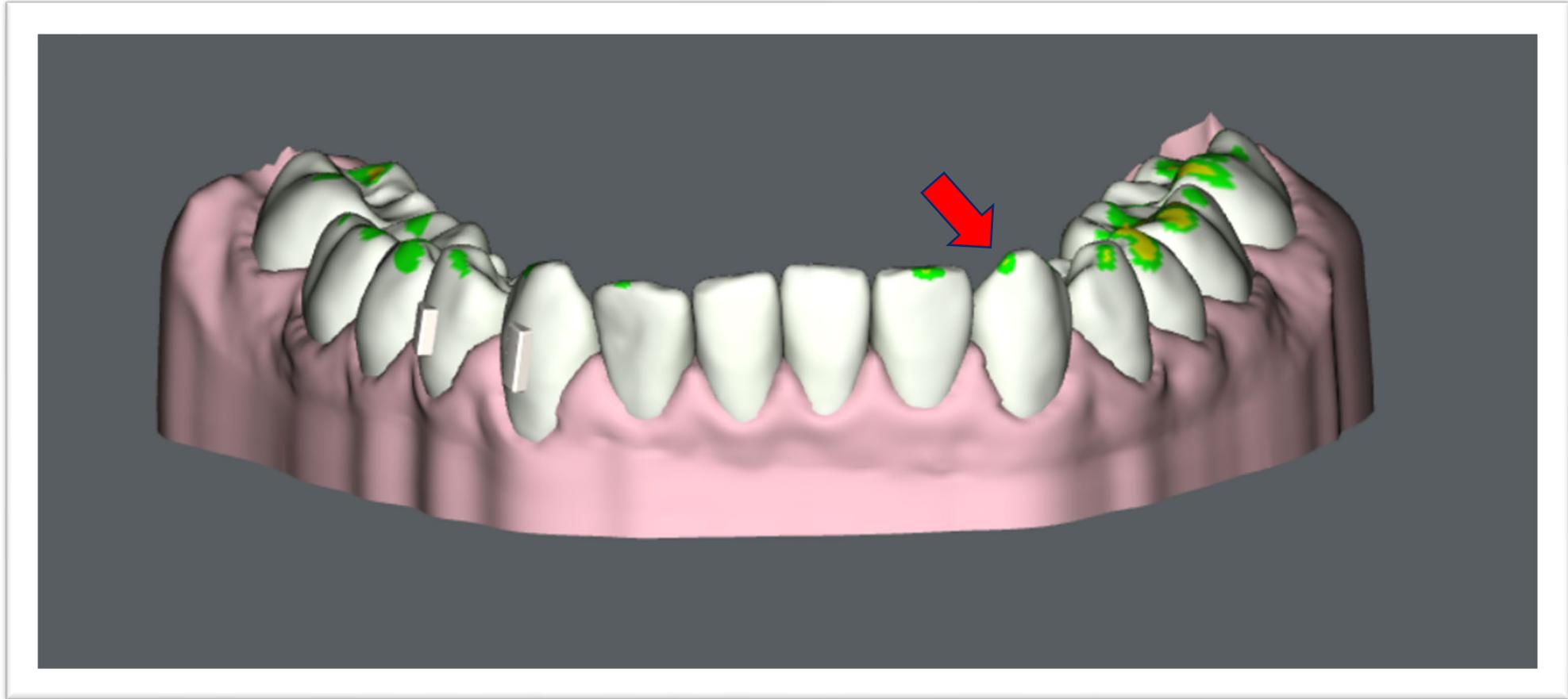

Notice light occlusal marking

Review SHO

Initial Stage

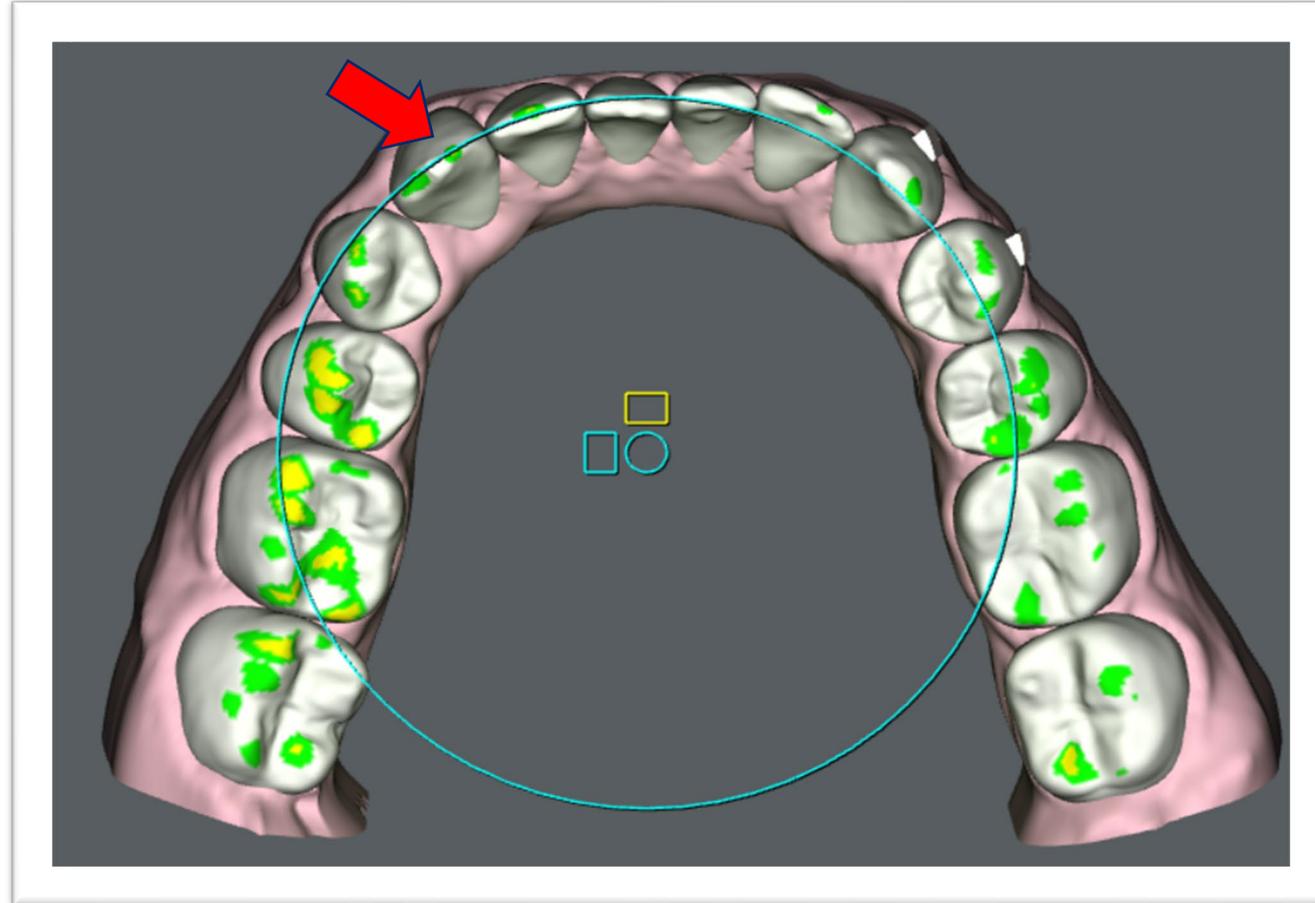

Notice light occlusal marking

Select Final Staged Model

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is "Jean Frost (R2N6)" and the current stage is "Setup & Staging 1 (Ready for Review)". A checklist on the left side shows the following items:


- Overview
- Setup and Staging Review
 - Labial View
 - Right Buccal View
 - Left Buccal View
 - Upper Occlusal
 - Lower Occlusal
 - IPR / Attachment / Trimline

The central 3D model shows a patient's teeth in a patient view. At the bottom, a navigation bar includes a timeline with markers at 1, 5, 10, 15, 20, and 21. The marker for 21 is highlighted with a red box, and a red arrow points to it from below. To the right of the timeline are icons for "Legend" and "Treatment Details".



Review SHO


Final Stage

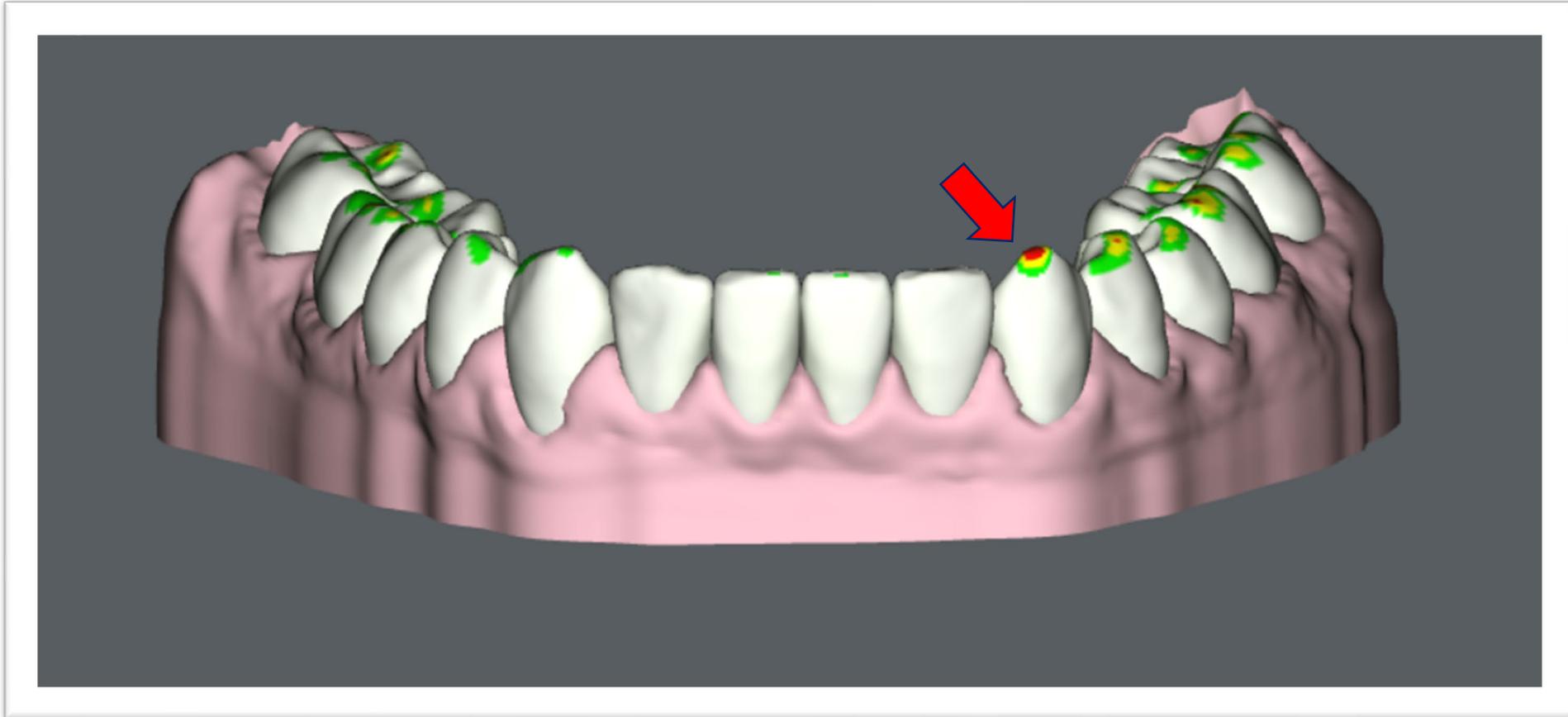

Notice purposeful over-corrected occlusal marking

Review SHO

Final Stage

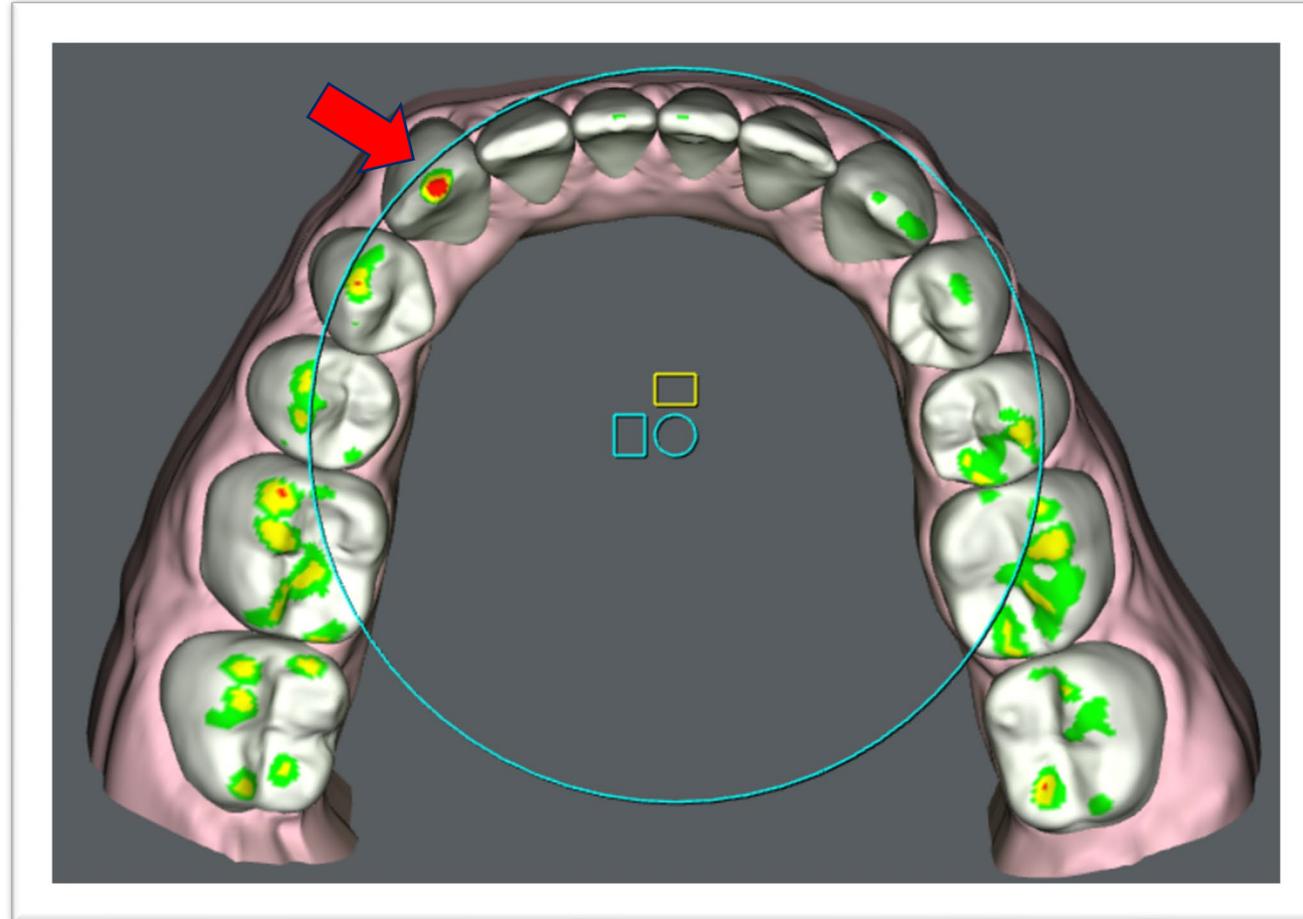

Notice purposeful over-corrected occlusal marking

Review THO

Transitory hyper-occlusion can occur when progressing through aligner treatment and will appear as red markings on the staged models. The patient may not be able to bite all the way down or may have sensitivity from prolonged hyper-occlusion like a high filling. Advise the patient at delivery this may occur and adjust if necessary.

Select Play Button

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is Jean Frost (R2N6) and the current stage is Setup & Staging 1 (Ready for Review). The main window shows a 3D model of a patient's teeth in a patient view. On the left, a 'Checklist and Overview' panel is visible, listing several review items, all of which are marked as complete with green checkmarks:

- Overview: Click the Notes Tab and review the notes from the Digital Lab (DL). Go to next step to begin reviewing the setup in different views, and review IPR
- Setup and Staging Review
 - Labial View
 - Right Buccal View
 - Left Buccal View
 - Upper Occlusal
 - Lower Occlusal
 - IPR / Attachment / Trimeline

Below the checklist, there is a 'Request Modifications' toggle switch and a button labeled 'Approve & Order Aligners'. At the bottom of the interface is a 'Camera (Patient View) Navigation' bar with various controls, including a play button. A red arrow points to this play button.

Review THO

Patient Database Tasks Denise Dodd (L5Y7) Setup & Staging 1 (To Review / Not ordered) SureSmile Aligner

Checklist and Overview


- Overview
- Setup and Staging Review
 - Labial View
 - Right Buccal View
 - Left Buccal View
 - Upper Occlusal
 - Lower Occlusal
 - IPR / Attachment / Trimline

Request Modifications

Approve & Order Aligners

Cancel Full Service Aligner Case

Camera (Object View) Navigation

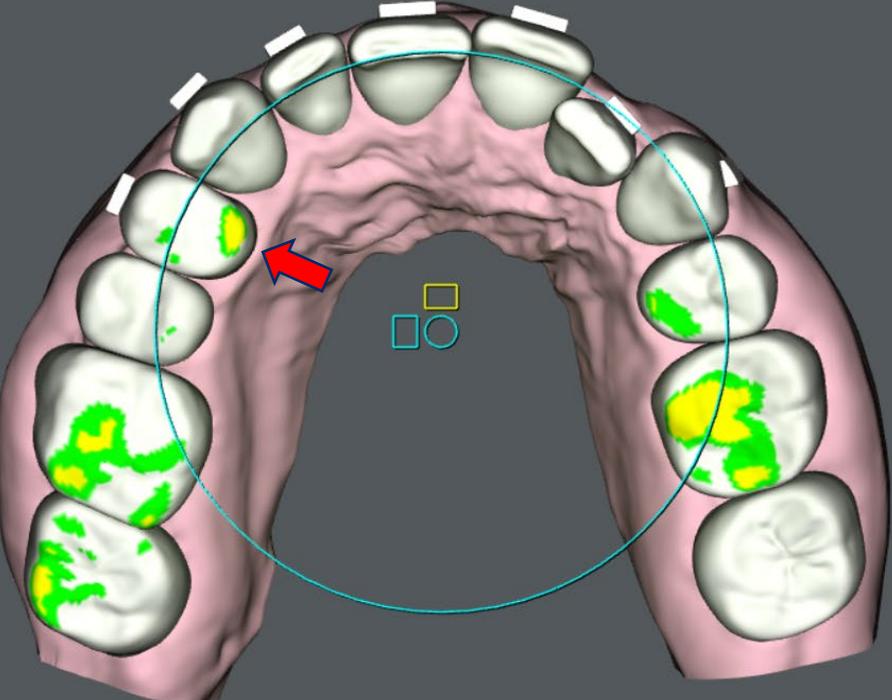

Legend Tooth movement

1 5 10 15 39



Review THO

The screenshot displays the SureSmile Aligner software interface. At the top, the patient name is Denise Dodd (L5Y7) and the case is Setup & Staging 1 (To Review / Not ordered). The main window shows a 3D model of the upper dental arch with occlusal markings in green and yellow. A red arrow points to the start of heavy occlusal marking on the left side. The left sidebar contains a Checklist and Overview with the following items: Overview (checked), Setup and Staging Review (checked), Labial View (checked), Right Buccal View (checked), Left Buccal View (checked), Upper Occlusal (checked), Lower Occlusal (checked), and IPR / Attachment / Trimline (not checked). Below the checklist is a 'Request Modifications' toggle and an 'Approve & Order Aligners' button. At the bottom, there is a 'Camera (Object View) Navigation' bar with a timeline from 1 to 39, a legend, and a 'Tooth movement' icon.



Notice start of heavy occlusal marking

Review THO

Patient Database Tasks Denise Dodd (L5Y7) Setup & Staging 1 (To Review / Not ordered) SureSmile Aligner


Checklist and Overview

- Overview
- Setup and Staging Review
 - Labial View
 - Right Buccal View
 - Left Buccal View
 - Upper Occlusal
 - Lower Occlusal
 - IPR / Attachment / Trimline

Request Modifications

Approve & Order Aligners

Cancel Full Service Aligner Case

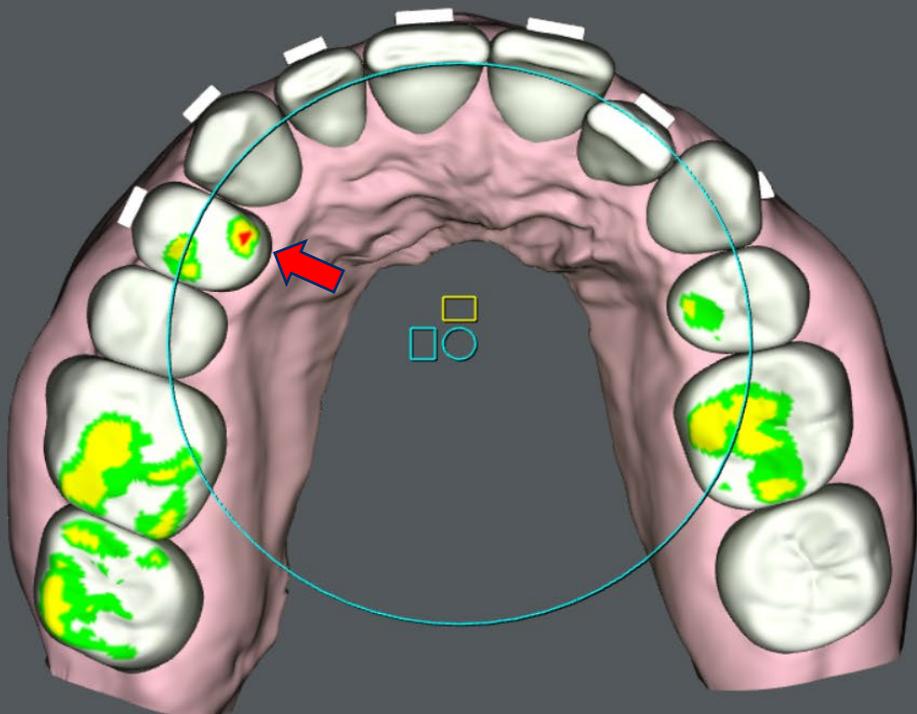

Camera (Object View) Navigation

1 5 10 12 15 39

Legend Tooth movement



Review THO

The screenshot displays the SureSmile Aligner software interface. At the top, the patient database shows "Denise Dodd (L5Y7)" and the case "Setup & Staging 1 (To Review / Not ordered)". The main view is a 3D model of the upper dental arch with occlusal markings in green and yellow. A red arrow points to the end of a heavy marking on tooth 18. On the left, a "Checklist and Overview" panel lists tasks: Overview (checked), Labial View (checked), Right Buccal View (checked), Left Buccal View (checked), Upper Occlusal (checked), Lower Occlusal (checked), and IPR / Attachment / Trimline (pending). Below the checklist is a "Request Modifications" toggle and an "Approve & Order Aligners" button. At the bottom, a "Camera (Object View) Navigation" bar shows a timeline with tooth numbers 1, 14, 15, 18, 20, 25, 27, and 39. A red arrow points to the 18 mark on this timeline. The interface also includes a "Legend" and "Tooth movement" section.

Notice end of heavy occlusal marking



CLEAR BLUE SMILES™

Step 3: Order Aligners

Step 3: Order Aligners

Follow complete aligner ordering workflow. Aligners should arrive within 5-7 business days. Schedule patient delivery visit two weeks from order date.

1. Select Approve and Order Aligners

The screenshot displays the SureSmile Aligner software interface. At the top, the patient information is Jackie Haumschild (H4J3) in the 'Setup & Staging 1' phase, which is 'To Review / Not ordered'. The main window shows a 3D dental model of the upper and lower teeth. A red circle highlights the central incisors, with a yellow square, a blue square, and a red circle overlaid on them. On the left, a 'Checklist and Overview' sidebar contains an 'Overview' section with instructions to review notes from the Digital Lab (DL) and a 'Setup and Staging Review' section with options for Labial View, Right Buccal View, Left Buccal View, Upper Occlusal, and Lower Occlusal. A 'Request Modifications' toggle is present. At the bottom of the sidebar, the 'Approve & Order Aligners' button is highlighted with a red arrow. Below the sidebar, there is a 'Cancel Full Service Aligner Case' button and a 'Camera (Patient View) Navigation' section with playback controls and a timeline from 1 to 14. The bottom right corner features a 'Legend' and 'Tooth movement' icon.

2. Select Complete Option



Choose your pricing option ✕

Complete Option

479.00 USD per patient, unlimited refinements for 3 years following first order
30.00 USD per retainer
Full service aligner: 150.00 USD
Refinements: free for 3 years

Select Option

15.00 USD per aligner
30.00 USD per retainer
Full service aligner: 150.00 USD
Refinements: 110.00 USD plus 15.00 USD per aligner

SureSmile® Retainers can be ordered any time after placing your initial aligner order. A reminder to "Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.

Trim Type & Height Close Next

3. Change Trim Height and Type

Choose your pricing option ✕

Complete Option

479.00 USD per patient, unlimited refinements for 3 years following first order
30.00 USD per retainer
Full service aligner: 150.00 USD
Refinements: free for 3 years

Select Option

15.00 USD per aligner
30.00 USD per retainer
Full service aligner: 150.00 USD
Refinements: 110.00 USD plus 15.00 USD per aligner

SureSmile® Retainers can be ordered any time after placing your initial aligner order. A reminder to "Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.



Select Trim Type

Trim Types & Heights
for aligner order

Aligner Trim Type Straight Trim Height

U Straight 1.0 mm

Scalloped

Aligner Trim Type Straight Trim Height

L Straight 1.0 mm

Scalloped

Note: Aligner Trim Type can be changed from Scalloped to Straight, but cannot be changed from Straight to Scalloped.

OK Cancel

"Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.

Trim Type & Height Close Next

Select Trim Height

Trim Types & Heights
for aligner order

Aligner Trim Type Straight Trim Height

U Straight 1.0 mm

Scalloped

Aligner Trim Type Straight Trim Height


L Straight 1.0 mm

Scalloped

Note: Aligner Trim Type can be changed from Scalloped to Straight, but cannot be changed from Straight to Scalloped.

"Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.

Trim Type & Height

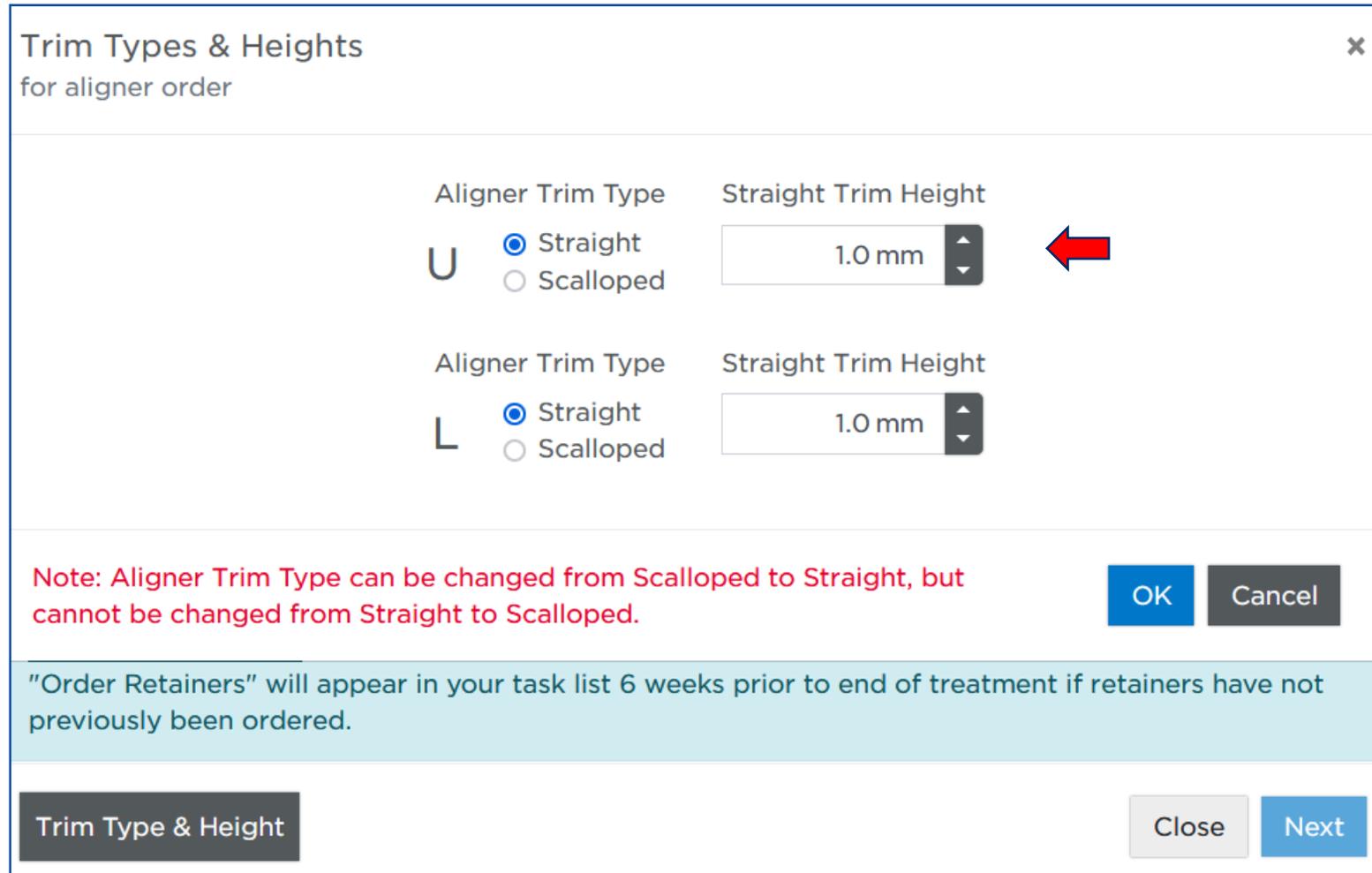

Select Trim Type

Trim Types & Heights
for aligner order

Aligner Trim Type Straight Trim Height

U Straight 1.0 mm

Scalloped

Aligner Trim Type Straight Trim Height

L Straight 1.0 mm

Scalloped

Note: Aligner Trim Type can be changed from Scalloped to Straight, but cannot be changed from Straight to Scalloped.

OK Cancel

"Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.

Trim Type & Height Close Next

Select Trim Height

Trim Types & Heights
for aligner order

Aligner Trim Type Straight Trim Height

U Straight 1.0 mm

Scalloped

Aligner Trim Type Straight Trim Height

L Straight 1.0 mm

Scalloped

Note: Aligner Trim Type can be changed from Scalloped to Straight, but cannot be changed from Straight to Scalloped.

"Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.

Trim Type & Height

OK Cancel

Close Next

4. Select OK

Trim Types & Heights
for aligner order

Aligner Trim Type Straight Trim Height

U Straight 1.0 mm

Scalloped

Aligner Trim Type Straight Trim Height

L Straight 1.0 mm

Scalloped

Note: Aligner Trim Type can be changed from Scalloped to Straight, but cannot be changed from Straight to Scalloped.

OK Cancel

"Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.

Trim Type & Height Close Next

5. Select Next

Choose your pricing option ✕

Complete Option

479.00 USD per patient, unlimited refinements for 3 years following first order
30.00 USD per retainer
Full service aligner: 150.00 USD
Refinements: free for 3 years


Select Option

15.00 USD per aligner
30.00 USD per retainer
Full service aligner: 150.00 USD
Refinements: 110.00 USD plus 15.00 USD per aligner

SureSmile® Retainers can be ordered any time after placing your initial aligner order. A reminder to "Order Retainers" will appear in your task list 6 weeks prior to end of treatment if retainers have not previously been ordered.

Trim Type & Height

Close **Next**

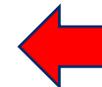


6. Select Place Order


Order Aligners ×

This order will incur a one-time charge of 629.00 USD including case (150.00 USD).

[Place Order](#)

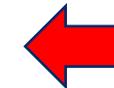

7. Select OK

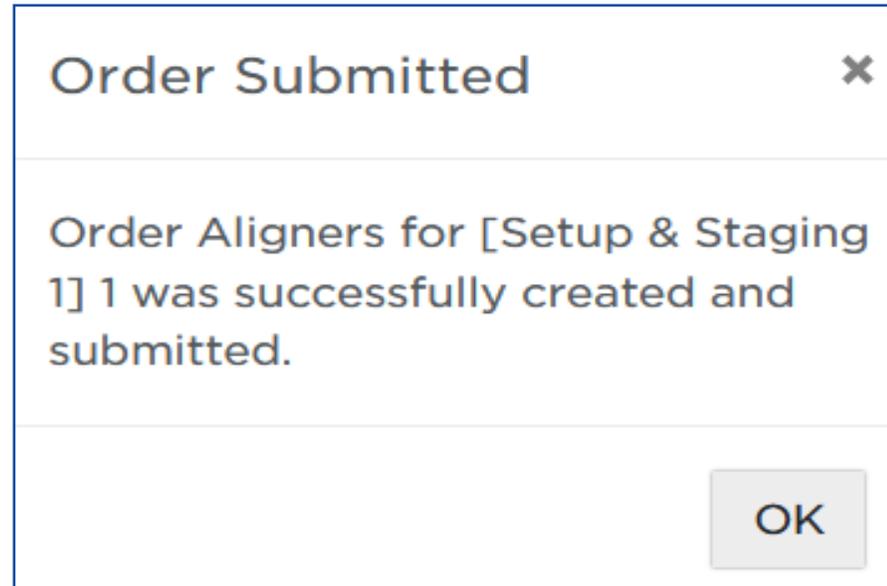

SureSmile Aligners



5-7 Business Days